

## ● How to Determine Glycan Profiles of Biopharmaceuticals from Peptide Mapping Data

A new and sensitive approach combining BioPharma Compass, PASEF and VIP-HESI

### Abstract

Glycosylation is a common critical quality attribute (CQA) of therapeutic proteins and needs to be characterized during product development. Typically, this analysis is conducted after enzymatic release and tagging of the glycans followed by fluorescence and MS detection. However, that approach loses information about the glycan position which is important for complex biologics such as fusion

proteins; and it requires dedicated method setups and experience far beyond the classical peptide mapping analysis.

Here we describe a workflow to identify glycan compositions directly from tryptic peptide maps [1] acquired with the timsTOF Pro and PASEF® in combination with analytical LC and the new VIP-HESI ion source. It employs a glycan search method from peptide mapping data and uses previously identified

aglycons as mass tags. Using this setup, we identified 36 different glycan compositions from NISTmAb tryptic digest which correlate well with those reported in the literature [2]. The results described here were generated directly from peptide mapping data without the need for a dedicated glycan laboratory workflow, thereby retaining the information about which peptide is glycosylated.

*Keywords:*  
PASEF, VIP-HESI,  
glycopeptides,  
peptide mapping,  
glycan composition,  
glycan site occupancy,  
timsTOF Pro

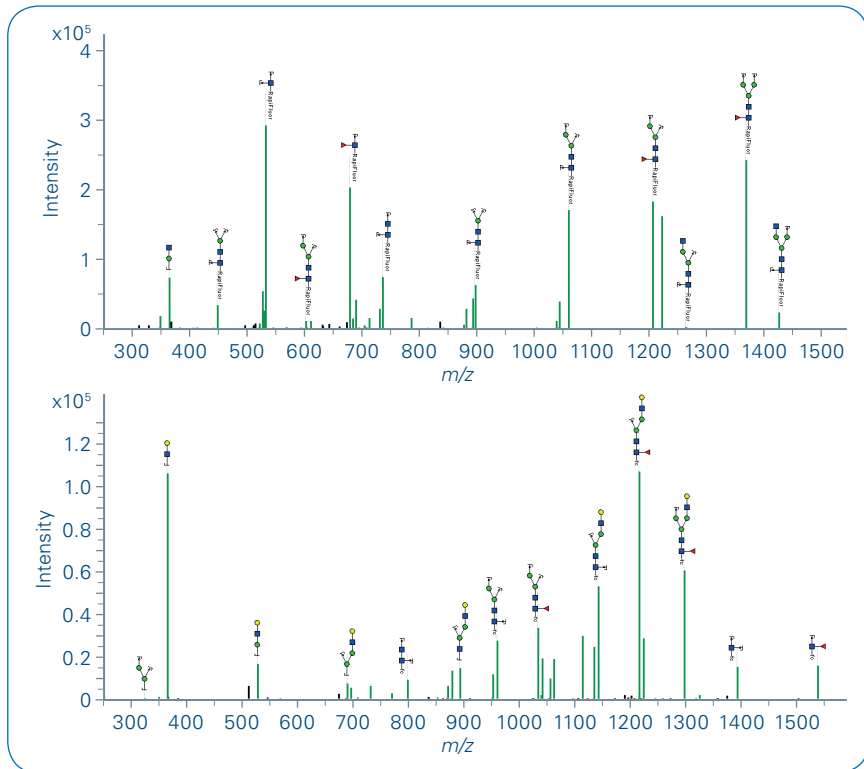


Figure 1. CID MS/MS spectra of G1F glycan from NISTmAb both dominated by  $\gamma$ -ions. Top: G1F with RapiFluor label (plus 311.175 Da) attached to the reducing end. Bottom: G1F with the peptide EEQYNSTYR (plus 1670.47 Da) assigned as "fc" attached to the reducing end.

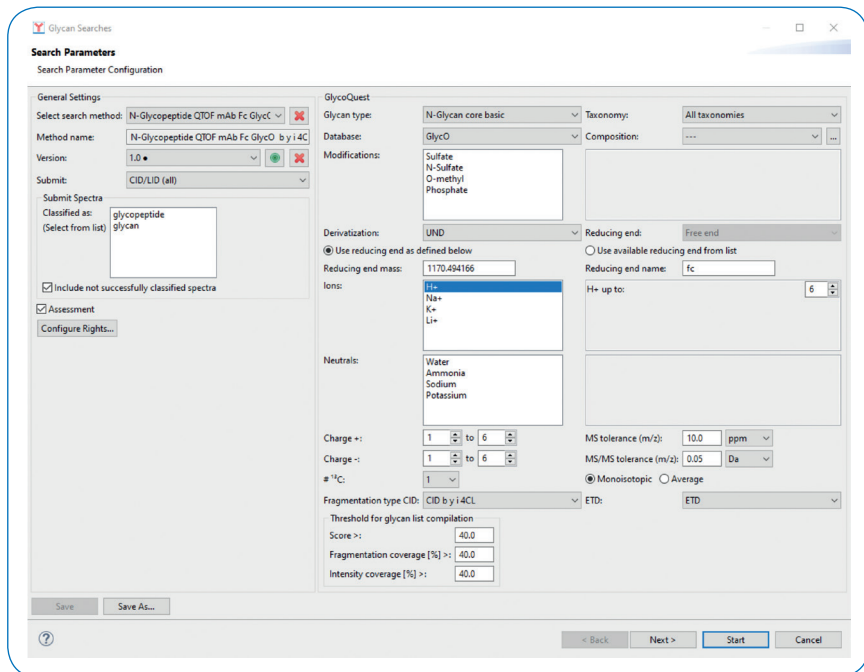


Figure 2. Glycan search method adapted for NISTmAb tryptic digest glycopeptide analysis. The monoisotopic reducing end mass is calculated from the tryptic glycopeptide sequence EEQYNSTYR. A new "Fragmentation type CID" is created, and high thresholds are used to reduce the number of false positives drastically listing only glycan compositions identified with high confidence. Especially the Intensity coverage value is important for the reliable identification of low abundant glycans.

## Introduction

Characterization of therapeutic proteins such as monoclonal antibodies requires a broad range of analyses including the full confirmation of the protein sequence as well as the detection and identification of post-translational modifications such as protein glycosylation. A common workflow to confirm the entire protein sequence is the peptide mapping approach, which combines proteolytic digests with RP-LC-MS/MS analysis in which, e.g., tryptic peptides are separated by reversed phase chromatography and further analyzed using high-resolution tandem mass spectrometry. N-linked glycans are subsequently analyzed in a second experiment comprising enzymatic glycan release, chemical labeling, LC separation on graphitized carbon or HILIC columns and MS analysis [2].

Glycopeptide analysis is typically applied in glycoproteomic experiments to identify the protein as well as the glycan using a specific fragmentation involving a collision energy stepping method [1]. The resulting glycopeptide MS/MS spectra provide both glycan and peptide fragments and allow to assign the peptide sequence as well as the corresponding glycan structure in one spectrum. Glycopeptides from therapeutic antibodies like NISTmAb, however, are derived from a tryptic peptide of known sequence and therefore can be assigned with reasonable confidence based on accurate mass alone.

When applied to glycopeptides, standard PASEF conditions used for peptide mapping preferentially cleave the glycosidic bond between the carbohydrate units, and subsequently peptide fragments are not observed. Such spectra are dominated by

y-type glycan fragments attached to the intact peptide and lower molecular weight b-type glycan fragments (without peptide moiety – i.e., the aglycon - attached). Similar fragmentation patterns are observed in CID spectra from released glycans carrying a fluorescence label like RapiFluor at the reducing end (Figure 1). In the approach reported here, the aglycon of a glycopeptide is treated in a glycan database search like a mass tag at the reducing end. The peptide tag serves the same purpose for the ionization of the glycan as fluorescent labels like RapiFluor in the analysis of labelled glycans – the peptide has an even stronger ionization propensity providing a high analysis sensitivity.

The timsTOF Pro with PASEF [3] further increases the capability to acquire high quality MS/MS spectra from glycopeptides with high dynamic range despite the almost non-existent reverse phase separation of the glycopeptides with identical aglycon. In this work we evaluated the identification characteristics of the glycopeptide approach on the timsTOF Pro linked to the VIP-HESI ion source (Vacuum Insulated Probe – Heated ESI) – a heated ion source optimized for applications requiring high sensitivity.

## Experimental

### Sample Preparation

NIST Monoclonal Antibody Material 8671 (Merck) was reduced using DTT and TFE (25%) for 150 min at 56 °C and alkylated using IAA. Trypsin (Promega) digestion was performed overnight.

### LC-MS Data Acquisition

Twenty  $\mu\text{g}$  tryptic digest were separated in a 60 min gradient on an

Acquity CSH C18 2.1 x 100 mm 1.7  $\mu\text{m}$  column (Waters) using an Elute UHPLC (Bruker). The UHPLC was interfaced with a timsTOF Pro mass spectrometer via a VIP-HESI ion

source and peptides were analyzed by PASEF using the standard proteomics acquisition method, adapted by lowering the precursor intensity threshold to 400 counts.

Row	Composition	Int.	Score	IntCov. [%]	FragCov. [%]	Rt [min]	m/z meas.	z	Mr calc.	Accession
1	Hex3HexNAc4dHex1-fc	161013	89.6	94	86	8.38	878.6905	3	2633.0386	2620
2	Hex4HexNAc4dHex1-fc	149800	85.3	96	76	8.30	1398.5562	2	2795.0914	13972
3	Hex4HexNAc3dHex1-fc	27556	68.5	90	52	8.30	1297.0142	2	2592.0121	13971
4	Hex3HexNAc3dHex1-fc	24036	83.4	92	76	8.30	1215.9866	2	2429.9592	14458
5	Hex5HexNAc2-fc	21720	50.8	46	56	7.97	802.6506	3	2404.9276	234
6	Hex5HexNAc4dHex1-fc	20358	81.9	92	73	8.22	986.7229	3	2957.1442	13970
7	Hex4HexNAc3NeuGc1dHex1-fc	19296	75.9	88	65	10.29	967.3748	3	2899.1024	38160
8	Hex6HexNAc4dHex1-fc	16997	72.3	85	61	8.22	1040.7406	3	3119.1971	13422
9	Hex3HexNAc2dHex1-fc	9804	74.4	91	61	8.30	1114.4468	2	2226.8799	1643
10	Hex3HexNAc3-fc	7498	67.9	72	64	8.22	762.3089	3	2283.9013	1865
11	Hex7HexNAc4dHex1-fc	7076	69.7	84	57	8.18	1094.7607	3	3281.2499	13423
12	Hex5HexNAc4NeuGc1dHex1-fc	6800	66.8	85	52	10.17	1089.0838	3	3264.2346	13922
13	Hex5HexNAc3dHex1-fc	6316	75.0	87	65	8.22	919.0319	3	2754.0649	14736
14	Hex5HexNAc3dHex1-fc	6316	78.1	92	66	8.22	919.0319	3	2754.0649	1248
15	Hex4HexNAc2-fc	5789	66.8	59	76	7.97	748.6333	3	2242.8748	2827
16	Hex5HexNAc4NeuGc2-fc	5564	54.3	73	40	10.05	1143.1029	3	3426.2757	13325
17	Hex4HexNAc4NeuGc1dHex1-fc	4869	78.4	86	72	10.25	1035.0672	3	3102.1818	13924
18	Hex3HexNAc5dHex1-fc	4835	64.6	57	73	8.46	946.3814	3	2836.1180	1549
19	Hex4HexNAc5dHex1-fc	4667	52.1	49	56	8.46	1000.3991	3	2998.1708	13904
20	Hex6HexNAc5dHex1-fc	3290	56.3	55	58	8.38	1108.4347	3	3322.2764	13909
21	Hex5HexNAc5dHex1-fc	3283	67.6	77	59	8.42	1054.4152	3	3160.2236	13392
22	Hex6HexNAc4NeuGc1dHex1-fc	2963	62.9	75	53	10.01	857.5796	4	3426.2874	13941
23	Hex6HexNAc3-fc	2484	71.5	84	61	8.01	924.3597	3	2770.0598	1710
24	Hex6HexNAc3NeuGc1-fc	2235	67.8	82	56	10.05	1026.7246	3	3077.1501	38059
25	Hex3HexNAc4-fc	2043	64.5	55	76	8.26	1244.4968	2	2486.9807	1857
26	Hex4HexNAc3-fc	1998	64.5	68	61	8.18	1223.9848	2	2445.9541	13944
27	Hex6HexNAc3NeuGc1dHex1-fc	1512	56.8	70	46	10.05	1075.4088	3	3223.2080	38101
28	Hex6HexNAc3dHex1-fc	1496	54.3	73	40	8.18	973.0463	3	2916.1177	1247
29	Hex6HexNAc2-fc	1280	69.2	77	62	7.81	856.6681	3	2566.9804	6296
30	Hex3HexNAc2-fc	1185	70.4	88	57	8.22	1041.4188	2	2080.8219	773
31	Hex5HexNAc3NeuGc1-fc	1133	65.6	75	58	10.17	972.7046	3	2915.0973	38094
32	Hex5HexNAc4dHex2-fc	1085	43.9	45	43	8.46	1035.4107	3	3103.2022	10137
33	Hex6HexNAc4dHex2-fc	1062	45.0	48	43	8.42	1089.4280	3	3265.2550	14487
34	Hex4HexNAc4dHex2-fc	971	63.7	67	61	8.54	981.3913	3	2941.1493	13206
35	Hex4HexNAc3NeuGc1-fc	743	59.4	76	46	10.29	918.6874	3	2753.0445	38143
36	Hex5HexNAc3-fc	574	51.4	48	55	8.14	1305.0067	2	2608.0070	13711

Figure 3. Result table of glycan compositions identified from NISTmAb tryptic glycopeptides sorted by absolute intensity (Int.). Hex=Hexose; HexNAc=N-Acetyl hexosamine; dHex=desoxyhexose; NeuGc= N-glycolylneuraminic acid; fc=reducing end mass of the tryptic Fc peptide EEQYNSTYR. Besides the Score from the glycan search software GlycoQuest, intensity coverage (IntCov. [%]) and fragmentation coverage (FragCov. [%]) are listed showing the high confidence level of the search results.

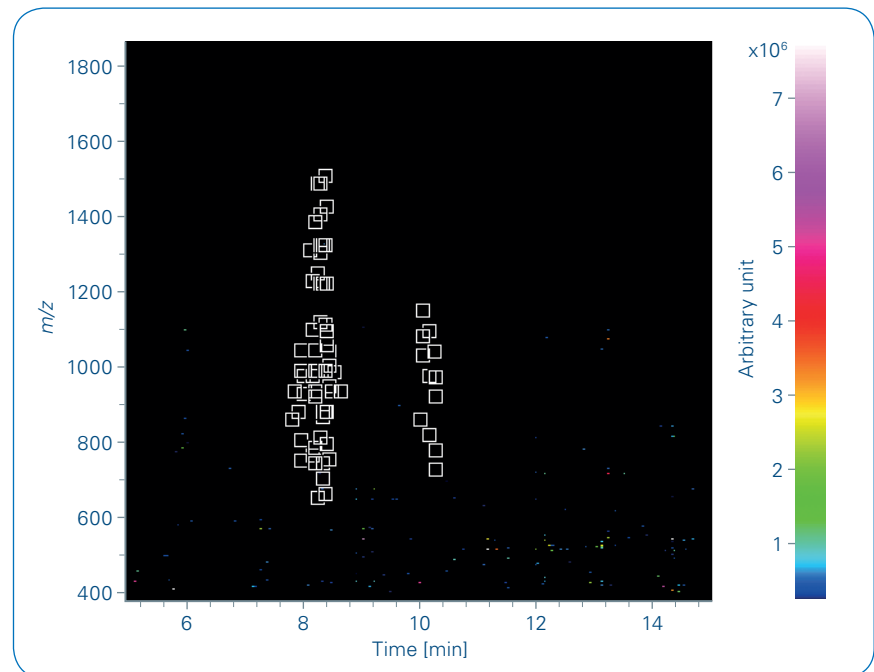


Figure 4. Survey view of identified glycopeptides at two different retention time ranges.

## Data Analysis

The raw data were initially processed using BioPharma Compass® 2021b (Bruker) with the Peptide Mapping method [Tutorial NIST\\_mAb](#) with small modifications resulting in classical peptide maps with comprehensive sequence coverage (not shown). Typically, the Fc-glycopeptides are eluting in a narrow retention time range early in the gradient. To reduce the number of MS/MS spectra submitted to the glycan search and computing time, the Rt range for data computing was limited to the glycopeptide elution range (5-15 min, Figure 4).

In a second processing step, all MS/MS spectra in the Rt range were submitted to a glycan search using the GlycoQuest search engine included within BioPharma Compass.

The Fragmentation Type *CID byi4CI* was initially defined in "Admin Preferences": Protocols/Protein-Scape/Glycomics/GlycoFragmentationType. Only b-, y- and internal ions are selected. *Maximum cleavages* are set to 4, *max crosslinks* to 0.

The tutorial glycan search method *N-glycan QTOF CID RapiFluor* was adapted subsequently to the analysis of the previously known/established tryptic NISTmAb Fc glycopeptide EEQYNSTYR using the following parameters (Figure 2):

1. Use Glyco as database as it is a rather condensed database containing all relevant glycans expected on a humanized IgG1
2. Define the Reducing end: Reducing end mass: 1170.494166 (peptide mass of EEQYNSTYR -H<sub>2</sub>O was calculated using the Sequence Editor in BioPharma Compass) Reducing end name "fc" for Fc tryptic peptide

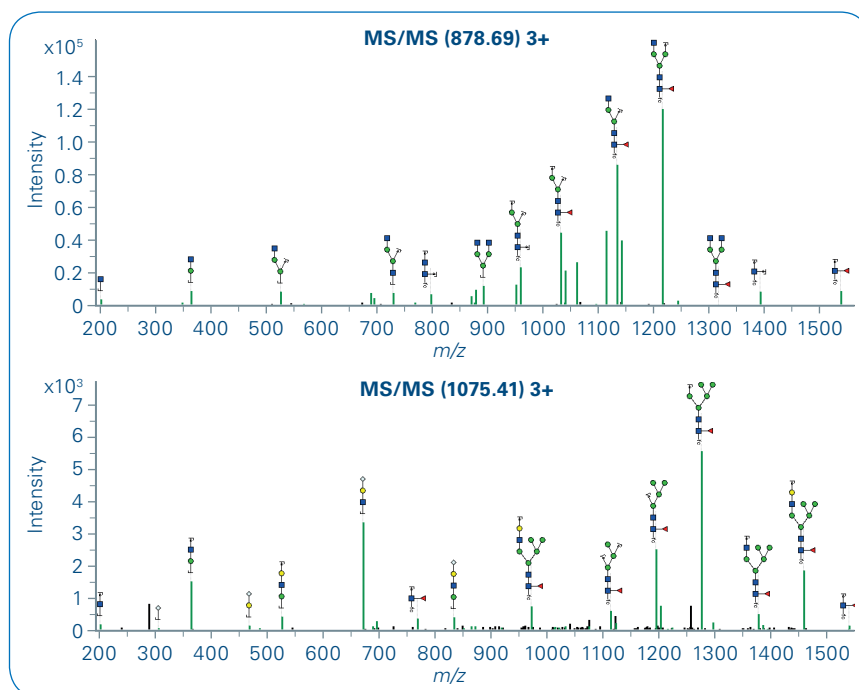


Figure 5. MS/MS spectra of two different glycan structures of different intensities with a dynamic range of 100 to 1. Top: G0F (h3n4f1; Hex3HexNAc4dHex1) attached to EEQYNSTYR with intensity of 161,013 a.u.. Bottom: Hybrid structure (h6n3f1g1, Hex6HexNAc3NeuGc1) with intensity 1,512 a.u..

3. Select Fragmentation Type *CID byi4CI*
4. Thresholds for result compilation were stringent to reduce the glycan list for confident search results: Score > 40, Fragmentation coverage [%] > 40, Intensity coverage [%] > 40

## Results

The glycan search of the MS/MS spectra yielded 36 specific glycan compositions of different glycan classes like complex, hybrid, and high mannose structures (Figure 3).

The glycopeptides were detected at 2 different Rt ranges: glycopeptides with neutral glycans and those with acidic glycan units (N-glycolylneuraminic acid) eluted after approx. 8 and 10 min, respectively (Figure 4).

The most intense glycan structure (G0F, h3n4f1) was detected with an absolute intensity of 161x10<sup>3</sup>, lower abundant compositions were

identified with intensities below 2x10<sup>3</sup> (Figure 5). These were mainly non-fucosylated neutral or acidic glycan structures or doubly fucosylated glycans like G1F2 or G2F2 (h4n4f2 or h5n4f2, Figure 6).

Compared to the results of a traditional analysis using released labeled glycans and HILIC separation [2], 27 out of 30 different compositions were found using our glycopeptide analysis method. The three missing compositions were h7n2, n7n3 and h8n5f1. The compositions h4n2, h4n3, h4n4f2, h5n3, h5n3g1, h5n4f2, h6n4g1 and h6n4f2 were additionally detected using the method described here (Figure 6). Differentiation of isomeric structures has not been further investigated.

The co-elution of glycopeptides always presents a challenge for LC-MS analysis, which typically results in much lower assignment rates of glycan compositions compared to labelled glycan analysis. Here, however, PASEF contributed

the required high acquisition speed (> 100 Hz MS/MS) and sensitivity (ion mobility time and space focusing), whilst sensitivity and signal intensities were given a further boost by the VIP-HESI ion source, which is heated to 400 °C. The combined effect was that 36 glycan compositions were identified from a single aglycon using mapping methods, which correlates well with data published by an expert group using a dedicated released glycan approach [2].

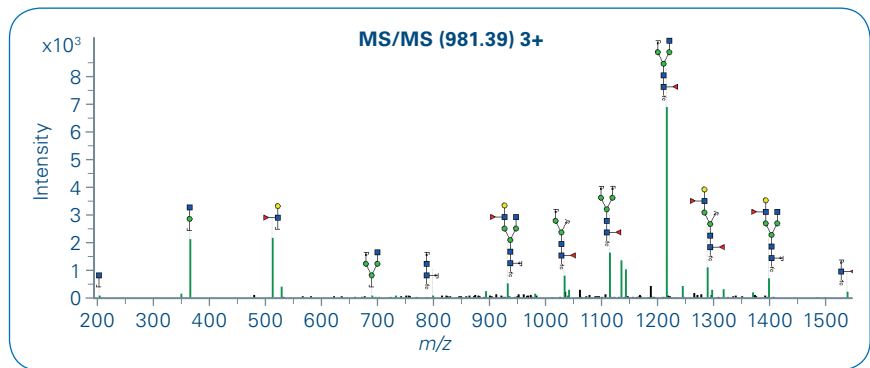


Figure 6. Additionally identified low abundant MS/MS spectrum of the di-fucosylated glycan composition h4n4f2 with intensity 971 a.u.. A score of 63.7 and intensity coverage of 67 underline the high stringency of the search result even for such a low abundant compound.

## Further reading

Bruker Application Note  
MT-132 (2021) Glycan Profiling and  
Sequencing SARS-CoV-2 RBD

Bruker Application Note  
1889046-VIP-HESI-2021-ebook

## Conclusion

- Thirty-six different glycan compositions were identified from a NISTmAb tryptic digest by dedicated glycoanalysis of a selected peptide using a conventional peptide mapping dataset and the GlycoQuest search engine in the BioPharma Compass software.
- The result compares favorably with previous studies using the established labelled-glycan approach [2]. However, a thorough isomer assessment still relies on MS/MS analysis and spectra library-based identification of labelled glycans.
- The PASEF technology was key to acquire high quality MS/MS spectra from coeluting glycopeptides covering a dynamic range of 100:1. The setup included analytical RP-LC separation and the new VIP-HESI ion source. This new approach enabled glycosylation analysis for non-experts.
- For more complex glycoproteins with multiple glycosylation sites, such as the SARS-CoV-2 S-glycoprotein [4], the approach can be repeated peptide-wise to obtain site-specific glycan composition profiles with great sensitivity and specificity.



## Learn More

You are looking for further Information?  
Check out the link or scan the QR code.

[www.bruker.com/timstofpro](http://www.bruker.com/timstofpro)



### References

- [1] Hinneburg H, Stavenhagen K, Schweiger-Hufnagel U, Pengelley S, Jabs W, Seeberger PH, Varón Silva D, Wührer M, Kolarich D (2016). *The art of destruction: Optimizing collision energies in Quadrupole – Time of Flight (QTOF) instruments for glycopeptide based glycoproteomics*. J Am Soc Mass Spectrom, DOI: 10.1007/s13361-015-1308-6
- [2] Hilliard M, Alley WR Jr, McManus CA, Yu YQ, Hallinan S, Gebler J, Rudd PM (2017). *Glycan characterization of the NIST RM monoclonal antibody using a total analytical solution: From sample preparation to data analysis*. mAbs, DOI: 10.1080/19420862.2017.1377381
- [3] Meier F, Brunner AD, Koch S, Koch H, Lubeck M, Krause M, Goedecke N, Decker J, Kosinski T, Park MA, Bache N, Hoerning O, Cox J, Räther O, Mann M (2018). *Online Parallel Accumulation–Serial Fragmentation (PASEF) with a Novel Trapped Ion Mobility Mass Spectrometer*. Mol Cell Proteomics **17**(12):2534-2545, doi.org/10.1074/mcp.TIR118.000900
- [4] Gstöttner C, Zhang T, Resemann A, Ruben S, Pengelley S, Suckau D, Welsink T, Wührer M, Dominguez-Vega E (2021). *Structural and Functional Characterization of SARS-CoV-2 RBD Domains Produced in Mammalian Cells*. Anal Chem **93**(17):6839-6847. doi: 10.1021/acs.analchem.1c00893. Epub 2021 Apr 19.

For Research Use Only. Not for Use in Clinical Diagnostic Procedures.

● **Bruker Daltonics GmbH & Co. KG**   **Bruker Scientific LLC**

Bremen · Germany  
Phone +49 (0)421-2205-0

Billerica, MA · USA  
Phone +1 (978) 663-3660

[ms.sales.bdal@bruker.com](mailto:ms.sales.bdal@bruker.com) – [www.bruker.com](http://www.bruker.com)