

## Errata Notice

This document contains references to PSS or Polymer Standards Service. Please note that PSS is now Agilent. This document will be republished as an Agilent document in the future.



A part of Agilent

# # 10277 - Column Application Note Characterisation of Glycerides

Modified biopolymer for production of emulsions in food-, cosmetic- and pharma industries. This GPC/SEC method is described in Ph. Eur. 2000.

## Experimental Setup

Mobile Phase:	Tetrahydrofuran
Stationary Phase:	PSS SDV
Flow rate [mL/min]:	1,00
Temperature [°C]:	25
Detection:	GPC1100 Refractive index
Calibration:	HPLC Calibration
Data processing:	PSS WinGPC Unity

## Recommendations for Sample Concentration


narrow PDI	
M 100 Da - 10 000 Da:	2 g/L
M 10 000 Da - 1 000 000 Da:	1-2 g/L
M > 1 000 000 Da:	0.5 g/L or less
broad PDI (>1.5)	
all molar masses:	3.0 - 5.0 g/L
Injection volume [ $\mu$ L]:	40




## Suitable Columns

low molecular weights:	S/N sda080505 and sda0860051e2
medium molecular weights:	-
high molecular weights:	-
ultrahigh molecular weights:	-

## Elugram overlay of triglyceride, diglyceride, and monoglyceride and glycerol

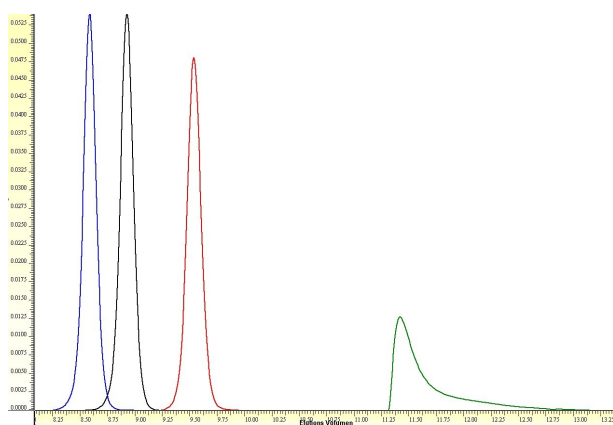

## Elugram of a technical mixture separation on PSS SDV

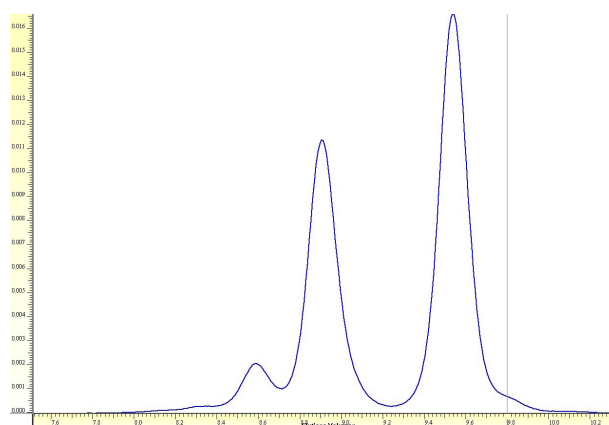

PSS Polymer Standards  
Service GmbH  
In der Dalheimer Wiese 5  
55120 Mainz | Germany

Phone +49 6131 96239-0  
Fax +49 6131 96239-11  
E-Mail [info@pss-polymer.com](mailto:info@pss-polymer.com)  
Web [www.pss-polymer.com](http://www.pss-polymer.com)

Polymer Standards  
Service-USA, Inc.  
160 Old Farm Rd, Suite A  
Amherst | MA 01002 | USA

Phone +1 413 835-0265  
Fax +1 413 835-0354  
E-Mail [pssusa@pss-polymer.com](mailto:pssusa@pss-polymer.com)  
Web [www.pss-polymer.com](http://www.pss-polymer.com)

RA45097.5465393519

5994-6257EN  
July 1, 2023