


Apple Juice with Added Standards

Rt- β DEXse™

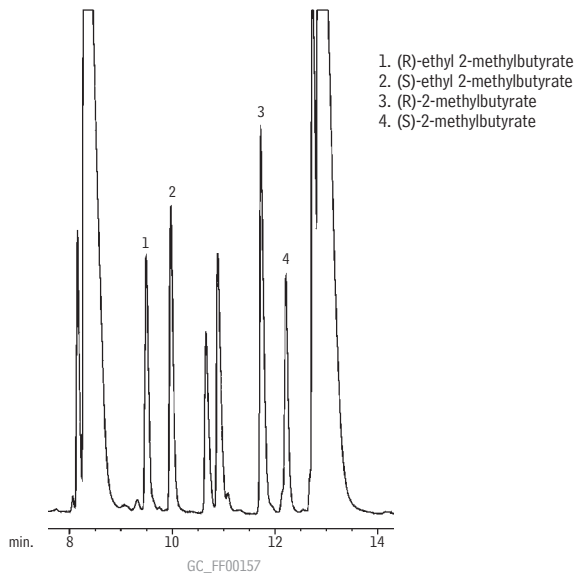

30m, 0.32mm ID, 0.25 μ m Rt- β DEXse™ (cat.# 13106)
1.0 μ L split injection
On-column conc. (standards): ~50ng each labeled enantiomer

Oven temp.: 40°C (hold 1 min.) to 220°C @ 2°C/min.
Inj. & det. temp.: 220°C
Carrier gas: hydrogen
Linear velocity: 80cm/sec.

Restek Corporation 110 Benner Circle Bellefonte, PA 16823
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • www.restek.com