


FAMES (Saw Palmetto)

Rtx[®]-Wax

Peak List	Conc. (mg/mL)
1. methyl caproate (C6:0)	0.4
2. methyl caprylate (C8:0)	0.4
3. methyl nonanoate (C9:0)	2.0
4. methyl caprate (C10:0)	0.4
5. methyl laurate (C12:0)	5.0
6. methyl myristate (C14:0)	2.0
7. methyl palmitate (C16:0)	2.0
8. methyl palmitoleate (C16:1)	0.4
9. methyl stearate (C18:0)	0.4
10. methyl oleate (C18:1)	5.0
11. methyl linoleate (C18:2)	1.0
12. methyl linolenate (C18:3)	0.4

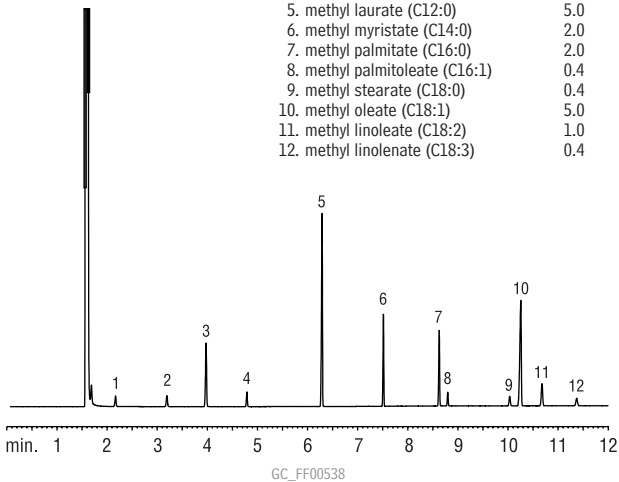

30m, 0.25mm, 0.25 μ m Rtx[®]-Wax (cat.# 12423)
1 μ L split injection of saw palmetto standard
Conc.: see peak list

Oven temp.: 120°C (hold 3 min.) to 220°C
at 20°C/min. (hold 12 min.)
Inj./det. temp.: 250°C/300°C
Carrier gas: helium
Linear velocity: 1mL/min. (34 cm/sec.)
Split ratio: 100:1