


FAMES (Ocean Nutrition Sample) on FAMEWAX

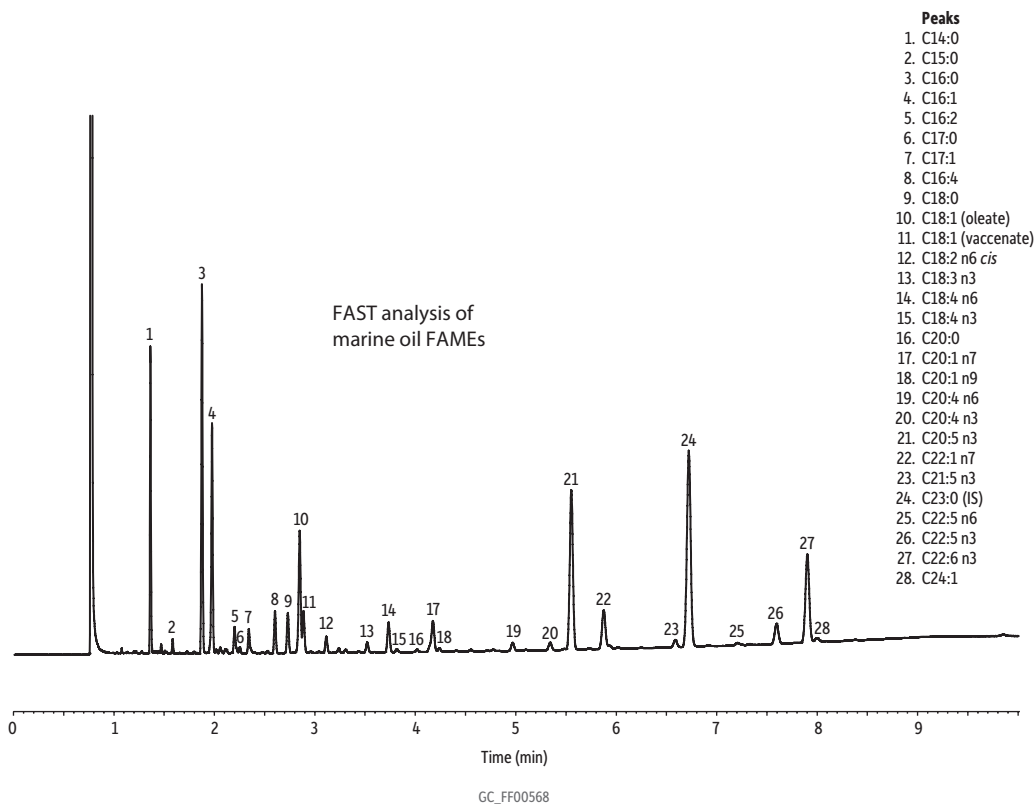

Column	FAMEWAX, 30 m, 0.32 mm ID, 0.25 µm (cat.# 12498)
Sample	Ocean nutrition sample
Conc.:	12 mg/mL total FAMES
Injection	
Inj. Vol.:	0.5 µL split (split ratio 150:1)
Liner:	3 mm split w/ wool
Inj. Temp.:	250 °C
Oven	
Oven Temp.:	195 °C to 240 °C at 5 °C/min (hold 1 min)
Carrier Gas	H ₂ constant flow
Linear Velocity:	62 cm/sec
Detector	FID @ 250 °C